


# PRECISION SCALE PRO

PRODUCT NUMBER 9432

Laboratory scales for medical and industrial applications

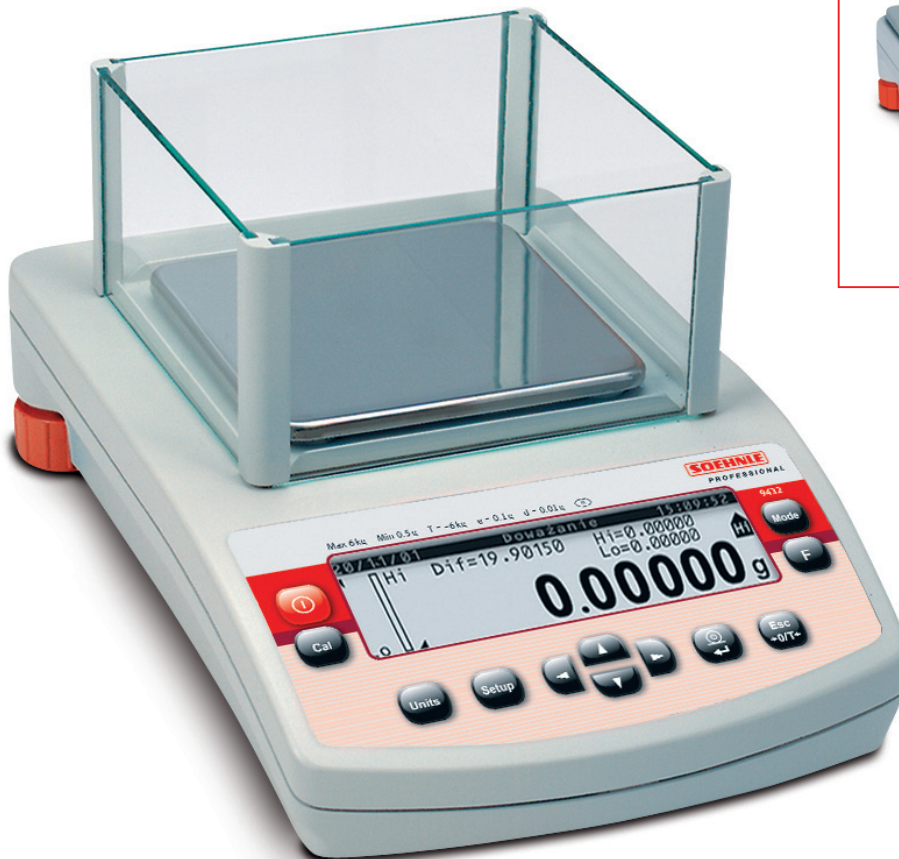

- > Weighing range: 200 g – 6,000 g
- > Division from 0.001 g
- > Large graphical and back-lighted display
- > Automatic internal calibration

# PRECISION SCALE PRO

PRODUCT NUMBER 9432

## Features



The series 9432 is the luxury class of Precision Scales and is perfect for fast definition of mass in laboratory- and industrial sectors. It has a weighing pan of stainless steel, a back-lighted graphical display and distinguishes itself by a high resolution. Data like Username, project name, date and time is programmable.

### Functions

- › Piece counting
- › Percentage weighing
- › Checking
- › Dosing
- › Recipes
- › Animal weighing
- › Density determination\*
- › Statistics
- › GLP procedures
- › Underfloor weighing
- › Programmable printouts

### Technical data

- › Operating temperature: + 15° till + 30° C
- › Power supply: Mains adapter 230 V AC, 9 V DC
- › Menu navigation: 6 languages selectable
- › Automatic internal calibration
- › Stabilisation time: 2 seconds
- › Interface RS 232
- › PS2 Interface to keyboard
- › Real time clock
- › Power supply included in delivery scope
- › Interface for an additional display possible

### Optional

- › Density determination set
- › Platform for underfloor weighing
- › Thermal printer with clock module
- › Thermal printer with statistics function
- › Tape printer
- › Wind deflector
- › Anti-vibration desk made of steel
- › PC software for transmitting and saving measuring data

## Ordering information

Order number	Weighing range	Division		Dimensions	Platform size
	g	g (d)	g (e)	mm	mm
9432.01.010*	200 / 2,000	0.001 / 0.01	0.01 / 0.1	335 x 210 x 80	128 x 128
9432.01.020*	750	0.001	0.01	335 x 210 x 80	128 x 128
9432.01.030	1,500	0.01	0.1	335 x 210 x 80	195 x 195
9432.01.040	4,500	0.01	0.1	335 x 210 x 80	195 x 195
9432.01.050	6,000	0.01	0.1	335 x 210 x 80	195 x 195

Your Soehnle Professional dealer:

[www.soehnle-professional.com](http://www.soehnle-professional.com)